

TRUECLOT BLEEDING CONTROL KIT SIMULATOR (1 GSW & 1 Laceration)


REV031824

TrueClot® Task Trainers are an excellent way to introduce your trainees and students to realistic hemorrhage control for **TCCC, TECC, TEMS, PHTLS, and general wound packing techniques for advanced hemostatic gauze training**

Dimensions Each:

- Measurements: H 17 in. x W 20 in. x D 7.5 in.
- Weight: 14 lb.

The Rescue Essentials Bleeding Control Instructors Kit contains everything you need to conduct a wet or dry wound packing training class. The TrueClot® Task Trainers are designed to simulate wound packing and demonstrate the use of advanced hemostatic gauzes with LUNA's TrueClot® Blood Simulant. The included, simulated large caliber gunshot wound, and knife/laceration Task Trainers can be packed with LUNA's Simulated Hemostatic Gauze while TrueClot® Blood Simulant is pumped into the wound using the attached 1-liter squeeze bottle. Realistic simulated blood clots form on the gauze at the point of bleeding demonstrating the function of hemostatic gauze.

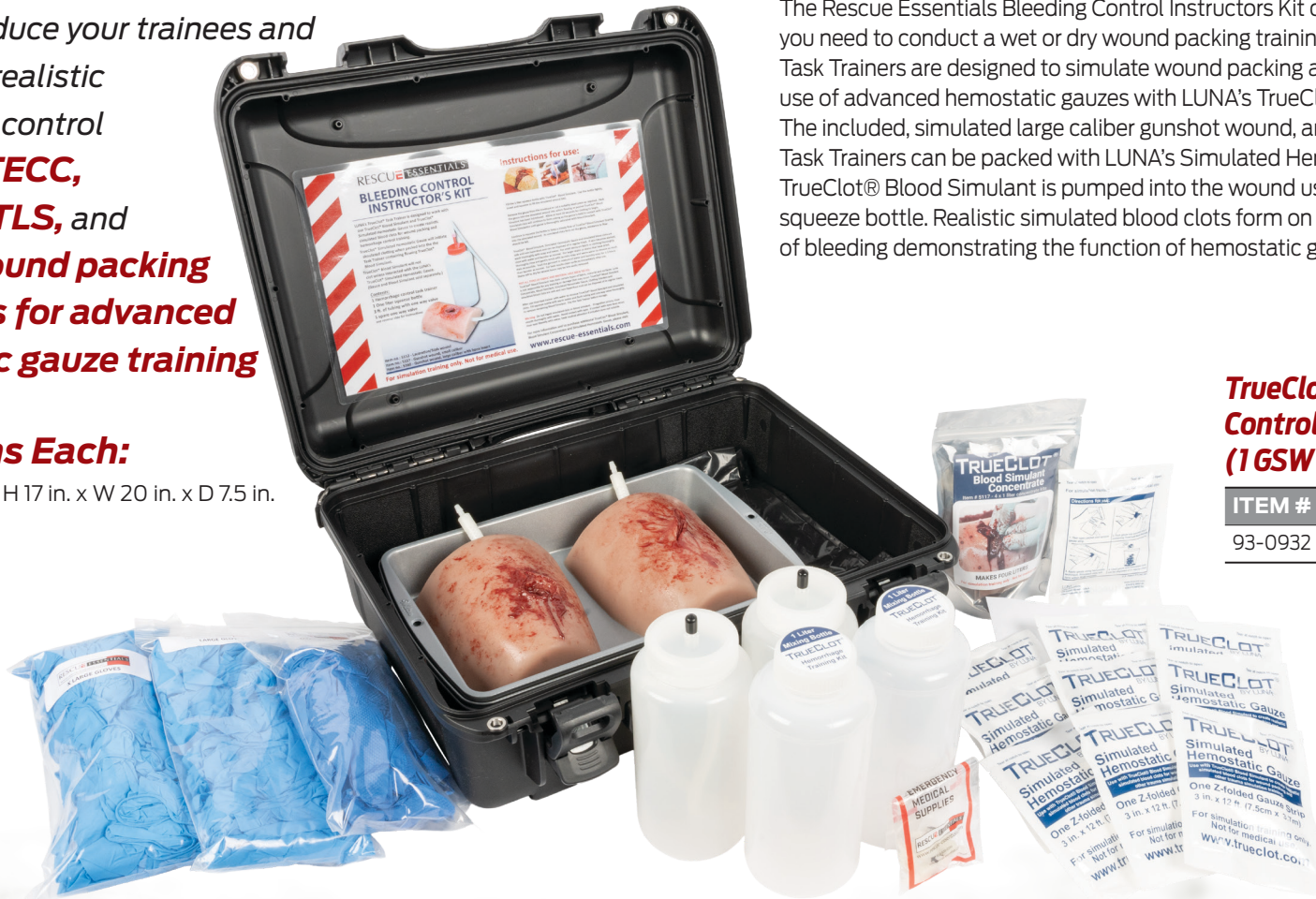

TrueClot® Bleeding Control Kit Simulator (1 GSW - 1 Laceration)

ITEM #	NSN #
93-0932	--



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A **MISSION®**

TRUECLOT BLEEDING CONTROL KIT SIMULATOR (1 GSW & 1 Laceration)

REV031824

Contents Include:

- 1 x TrueClot® Task Trainer with Bone and Large Caliber Wound (GSW)
- 1 x TrueClot® Task Trainer with Knife/Laceration Wound
- 2 x Bleeding Bottles with Tubing and Valves
- 2 x Blood Mixing & Storage Bottles
- 1 x TrueClot® Replacement Flow Valves, 5 pack
- 4 x TrueClot® Blood Simulant Concentrate Packs (makes 1L of Blood Simulant each for 4 liters total)
- 10 x TrueClot® Simulated Hemostatic Gauze, Z Fold, 12ft
- 2 x Wound Pans for Blood Containment
- 12 x Absorbent Underpad Chux
- 2 x Disposal Bags for Used Training Supplies
- 20 pr Nitrile gloves, Medium
- 20 pr Nitrile gloves, Large
- 20 pr Nitrile Gloves, X large
- 1 x Storage & Transport Case



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A **MISSION®**