

M-FAK LCL *Mini First Aid Kit*

NAR®'s most compact, versatile Mini First Aid Kit, now available in laser cut nylon!

North American Rescue's Mini First Aid Kit (M-FAK) - now featuring Laser Cut Laminate (LCL) MOLLE - was designed to be the most compact, multi-use IFAK for delivering immediate critical care for penetrating, blast or other traumatic injuries. At half the size of our popular T.O.R.K. platform, the M-FAK is currently the smallest platform available in our product offering. Despite its reduced size, this kit comes fully loaded with the critical medical equipment most requested by First Responders operating in the line of duty. Treatment pack out options include Basic, Basic with Combat Gauze® LE and Advanced with Combat Gauze® LE versions.

This Individual First Aid Kit's super compact, rugged nylon bag opens via 550 cord zipper pulls from the side in a clam shell configuration utilizing (2) two main sleeves that open on both ends for easy access and features multiple elastic loops for secure gear organization. The bag's exterior includes a MOLLE backing for mounting in the standard vertical position as well as a 3 inch internal sleeve for accommodating horizontal attachment to your duty belt.

The vertical mount can be set to open left-to-right or right-to-left based on shooter preference, while the horizontal mount on a belt allows opening directly to your C-A-T® tourniquet. Also included are a loop patch on front of the bag for custom user labeling and (3) three MOLLE strips on the front, and (2) two MOLLE strips on the spine's edge to attach additional gear such as TQ/holders, tactical lights, or knives (sold separately) if needed.



NORTH AMERICAN RESCUE®

www.NARescue.com • 888.689.6277

M-FAK LCL

COLOR	CONFIGURATION	ITEM #
BLK	Basic	80-1667
COY	Basic	80-1668
OCP	Basic	80-1669
ODG	Basic	80-1670
RED	Basic	80-1671
RGR	Basic	80-1672
BLK	with Combat Gauze® LE	80-1673
COY	with Combat Gauze® LE	80-1674
OCP	with Combat Gauze® LE	80-1675
ODG	with Combat Gauze® LE	80-1676
RED	with Combat Gauze® LE	80-1677
RGR	with Combat Gauze® LE	80-1678
BLK	Adv. w/ Combat Gauze® LE	80-1679
COY	Adv. w/ Combat Gauze® LE	80-1680
OCP	Adv. w/ Combat Gauze® LE	80-1681
ODG	Adv. w/ Combat Gauze® LE	80-1682
RED	Adv. w/ Combat Gauze® LE	80-1683
RGR	Adv. w/ Combat Gauze® LE	80-1684

Features:

- NAR's most compact, versatile Mini First Aid Kit - now with laser cut laminated MOLLE
- Contains the First Responder's most requested critical point-of-wounding medical equipment for treating penetrating, blast or other traumatic injuries in the line of duty
- The LCL IFAKS are double spaced for the MOLLE attachment points on the rear of each IFAK for maximum adjustability and super secure mounting to any platform
- Half the size of our popular T.O.R.K. platform

REV072424

Nylon Colors:

BLK

COY

OCP

ODG

RED

RGR

PRODUCTS
WITH A **MISSION®**

M-FAK LCL *Mini First Aid Kit*

REV072424

Content options:

BASIC




QTY 1 C-A-T® Tourniquet
QTY 1 ETD® 4 in. Dressing
QTY 1 NAR® Gauze, Wound Packing
QTY 1 HyFin® Vent Compact, Twin Pk
PAIR 1 Responder Black Gloves, LG

with COMBAT GAUZE®



QTY 1 C-A-T® Tourniquet
QTY 1 ETD® 4 in. Dressing
QTY 1 Combat Gauze® LE
QTY 1 HyFin® Vent Compact, Twin Pk
PAIR 1 Responder Black Gloves, LG

ADVANCED

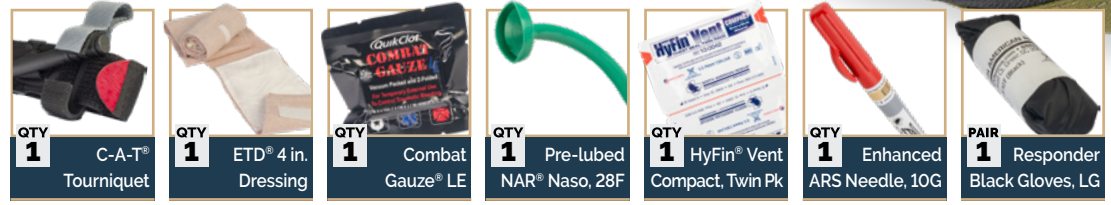

QTY 1 C-A-T® Tourniquet
QTY 1 ETD® 4 in. Dressing
QTY 1 Combat Gauze® LE
QTY 1 Pre-lubed NAR® Naso, 28F
QTY 1 HyFin® Vent Compact, Twin Pk
QTY 1 Enhanced ARS Needle, 10G
PAIR 1 Responder Black Gloves, LG

Specifications:

- L 6 in. x W 3 in. x D 3 in.
- Weight: 13 oz



NORTH AMERICAN RESCUE®
www.NARescue.com • 888.689.6277

PRODUCTS
 WITH A **MISSION®**